



Sun Gait Trainer

Height Adjustable Frame

Multi-adjustable Accessories

Light-weighted Sturdy Material



The Sun Gait Trainer is a mobility device crafted to provide support and assistance to youngsters who face challenges with their mobility.

Beneland Comfort LTD. / info@beneland.cn / <https://beneland.net>



Feature

The Sun Gait Trainer offers sturdy support, adjustable height, and customizable features to meet the unique needs of each user across various developmental stages.



One Frame, Multiple Possibilities

With a range of accessory options, the Sun Gait Trainer offers flexible anterior and posterior configurations for different stages of gait training.



Posterior Use



Anterior Use





Height Adjustable Frame



S: 5 height levels
M: 6 height levels
L: 8 height levels



Customization

The frame of the Sun Gait Trainer can be customized with additional accessories such as handles, backrests, and headrests according to the user's specific needs and circumstances. The Sun Gait Trainer features adjustable height, allowing it to adapt and grow with your child for continuous support and comfort.

User Height





Accessories




Headrest

The headrest provides stable and gentle support to the user's head, allowing the user to maintain correct and comfortable position while walking. Also, the headrest is crafted with antibacterial and breathable materials, providing a hygienic environment and preventing moisture buildup.

The adjustable features ensure a customized fit for each user.



* The headrest needs to be used with the backrest. It is not able to be attached to the frame directly.

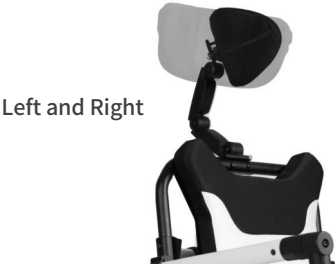

Left and Right



Up and Down



Front and Back

Backrest



The backrest provides support to the user's trunk during gait training, preventing backward falls. With three adjustable backrest angles (80°, 85°, and 90°), a gentle forward tilt helps move the user's center of gravity forward, supporting a smoother start to gait training.

Handles

Swivel Armrest with Handlegrip



The handledrip is fixed with the armrest, ensuring the handle rotates together with the armrest for better control and support during gait training.



Swivel Armrest with Adjustable Handlegrip



Swivel forearm support with a 360° rotating, angle-adjustable handlegrip for flexible positioning and improved comfort during use.



Upper Trunk Support Pad



The upper trunk support pads stabilize the user's upper body and provide sufficient support.

When used together with the backrest and saddle seat, no additional straps are required for fixation.



* The upper trunk support pad needs to be used with the backrest. It is not able to be attached to the frame directly.

Wrap-Around Trunk Support



The Wrap-Around Trunk Support stabilizes the user's trunk and provides reliable support. Its height and tilt angle are adjustable, and it can be used in both anterior and posterior modes of the Sun Gait Trainer.



Lower Trunk Support Strap



The Lower Trunk Support Strap stabilizes the user's core and helps prevent forward leaning or falling during training. It has a maximum weight capacity of 75 kg and is intended for users with mild symptoms and sufficient trunk control.



Leg Separator



The Leg Separator prevents the feet and lower legs from crossing the midline, guiding the child towards a more neutral and natural walking pattern.



Ankle Prompts



Ankle Prompts are designed to prevent the crossing of legs or scissoring during walking. They attach to the bottom of the gait trainer frame and can be adjusted in length for customized foot separation.



Posterior Saddle Seat



* The posterior saddle seat needs to be used with the backrest. It is not able to be attached to the frame directly.

Designed for posterior gait training, the adjustable and foldable saddle seat offers support while training and a convenient resting position for children when fatigued. Its height can be adjusted to suit different users, providing a flexible and comfortable option for breaks.

Anterior Support Seat



Designed for anterior gait training, the Anterior Support Seat helps support the user during transfers and setup. It allows caregivers to position the child securely before fitting upper body supports and other accessories. During training, it can also provide short-term resting support, but is not intended for full weight-bearing seating.



Multifunctional Castors

The Multifunctional Castors are designed with anti-slip materials to ensure safety and stability during use.



Front Castors

Equipped with swivel locks that can be set to allow forward rolling only, enhancing control and direction.



Rear Castors

Feature swivel locks for forward rolling, along with braking and speed control functionalities for added safety and customization.

For a straight walking path, lock all four swivel locks.

For stability with maneuverability, lock only the rear swivel locks to allow for controlled turning.



U Style Harness



The U style harness provides support to user's upper body, helping the user maintain an upright posture and improving upper body control during walking and gait training.

By securing the upper body area, the harness enhances the user's stability and balance, reducing the risk of falls and improving overall walking confidence.



Hip Positioner




The hip positioner provides support to the pelvis and hip area, helping the user maintain proper alignment and stability during walking and gait training exercises.

The hip positioner is designed to fit comfortably and securely, providing the necessary support without impeding the user's movement.

* The U style harness and hip positioner / pelvic support belt need to be used with the backrest. They are not able to be attached to the frame directly.




Customize and combine accessories freely to perfectly cater to each child's unique needs.

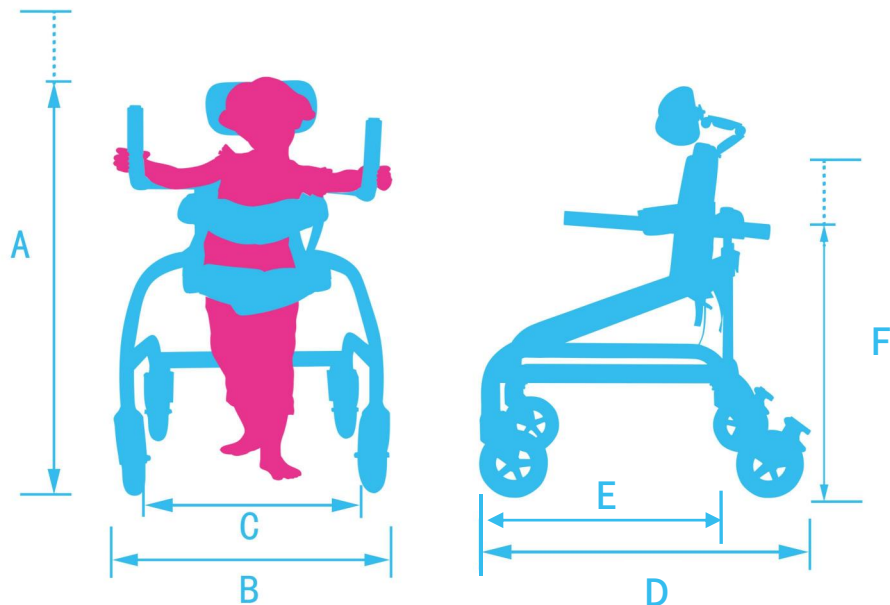


Item Dimensions	Small	Medium	Large
Overall Height Including The Headrest (A)	89 - 114 cm	104 - 139 cm	123 -164 cm
Overall Width (B)	55.5 cm	60 cm	62.5 cm
Inner Width (C)	44 cm	49 cm	52 cm
Overall Depth (D)	86 cm	90.5 cm	92 cm
Inner Depth (E)	62 cm	66 cm	70 cm
Height From Bottom to Handle (F)	53.5 - 67.5 cm	69 - 90 cm	86 - 113 cm
User Height	100 -125 cm	115 - 150 cm	134 - 175 cm

U Style Harness

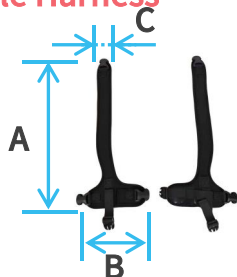

Item Dimensions	Small	Medium	Large
A	44.5	49 cm	54.5 cm
B	19 cm	21.5 cm	25 cm
C	5.5 cm	5.5 cm	7 cm

Hip Positioner

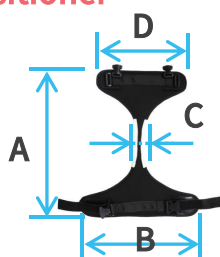

Item Dimensions	Small	Medium	Large
A	54 cm	59 cm	67 cm
B	40.5 cm	50.5 cm	58.5 cm
C	3.5 cm	4 cm	4 cm
D	30 cm	35.5 cm	48.5 cm
Waist	45.5 - 68.5 cm	66 - 81 cm	78.5 -101.5 cm